

DAVIES[®] GOLD

Dermatome and electric motor system

The Davies[®] Gold dermatome system is a precision engineered set with a proven track record for quality and dependability. The Duplex head provides multiple cutting widths and is offered with a variety of accessories. The head is designed for use with the CM2 controller and Micro motor. However, it can be driven by alternative motor systems. - Please ask for further details.



- Compact and lightweight
- Wide range of cutting depths
- Selection of cutting widths
- Robust, proven design
- Dermatome head, motor unit and cable are fully sterilisable
- Very accurate and consistent cutting performance
- Powerful, variable speed, high torque motor unit
- CE mark conformance
- Easily maintained and serviced

System specification

Dermotome specification

Cutting width (mm)	73
Cutting thickness range-7 settings	0.08 - 0.38 mm
Operating speed range (RPM)	2000 - 15000
Weight (Kg)	0.32

Optional

Reduction sole plates	Part No.
Duplex to 48 mm cutting width	8540
Duplex to 23 mm cutting width	8550

Micro motor and power supply

The dermatome head fits quickly onto the universal micro motor which is powered by the variable speed CM2 controller. The motor can be controlled either by the footswitch or directly from the controller itself. Please specify voltage and plug configuration when ordering.


Power supply specification

Input voltage:	230 (120) VAC/50 (60) Hz
Power:	130VA
Dimensions (mm):	100 x 85 x 165
Weight (Kgs):	1.68/0.12 controller/motor



Sterilisation Container

This purpose designed compatible delivery system will ensure your instruments are maintained in good condition at all times. The container is 100% sterilisable by steam, dry heat ('flash') to 160C, ethylene oxide (ETO) or gamma irradiation.

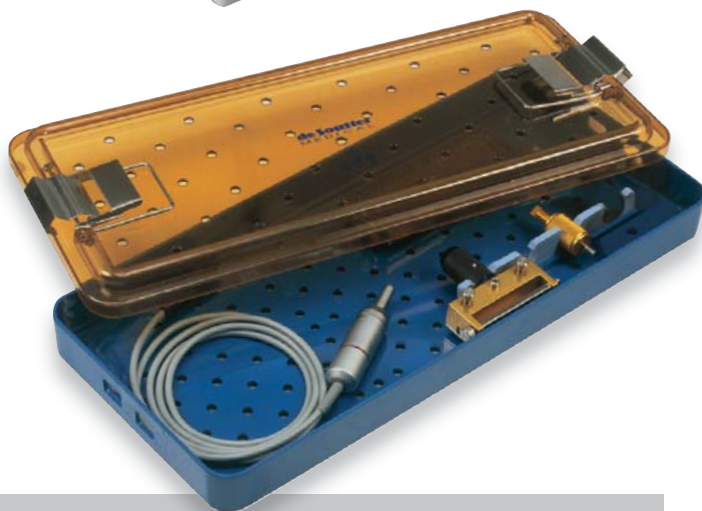

How to order

System components	230V - UK plug	230V - Euro plug	120V - US plug
GD-103 Duplex head complete	222994	222994	222994
CM2 motor controller and footswitch	1280884	1280994	1281004
Duplex blades (pack of 5)	S11-103	S11-103	S11-103
Sterilisation container	8570	8570	8570
Complete system ordering code	GDA-103/230UK	GDA-103/230EU	GDA-103/120US

De Soutter Medical Limited

Halton Brook Business Park, Weston Road, Aston Clinton
Aylesbury, Bucks HP22 5WF
tel: (0)1296 634 000 fax: (0)1296 634 033
email: info@de-soutter.com web: www.de-soutter.com

De Soutter Medical Australia Pty Ltd

2/12-14 Apollo Drive, Hallam, Victoria 3803
tel: 03 9702 4441 fax: 03 9702 4484
email: australia@de-soutter.com
web: www.de-soutter.com

deSoutter
MEDICAL