

# SPECIFICATION AND PART NUMBERS

## Irrigation pump kits



Irrigation tube fitted to a De Soutter Medical HST-100. For more information on this model, see catalogue no. 164. For information on Model EPV, see catalogue no. 170.



Description	Order code
Irrigation tube set - bur	T10-100
Irrigation tube set - HST-100	T10-102

## Wash Assist hose assemblies




Instrument type	Hose length and order codes		
	3m	4m	5m
EPV, HST	632333	632343	632353

## IPS-220 pump sensor



Order code 11570

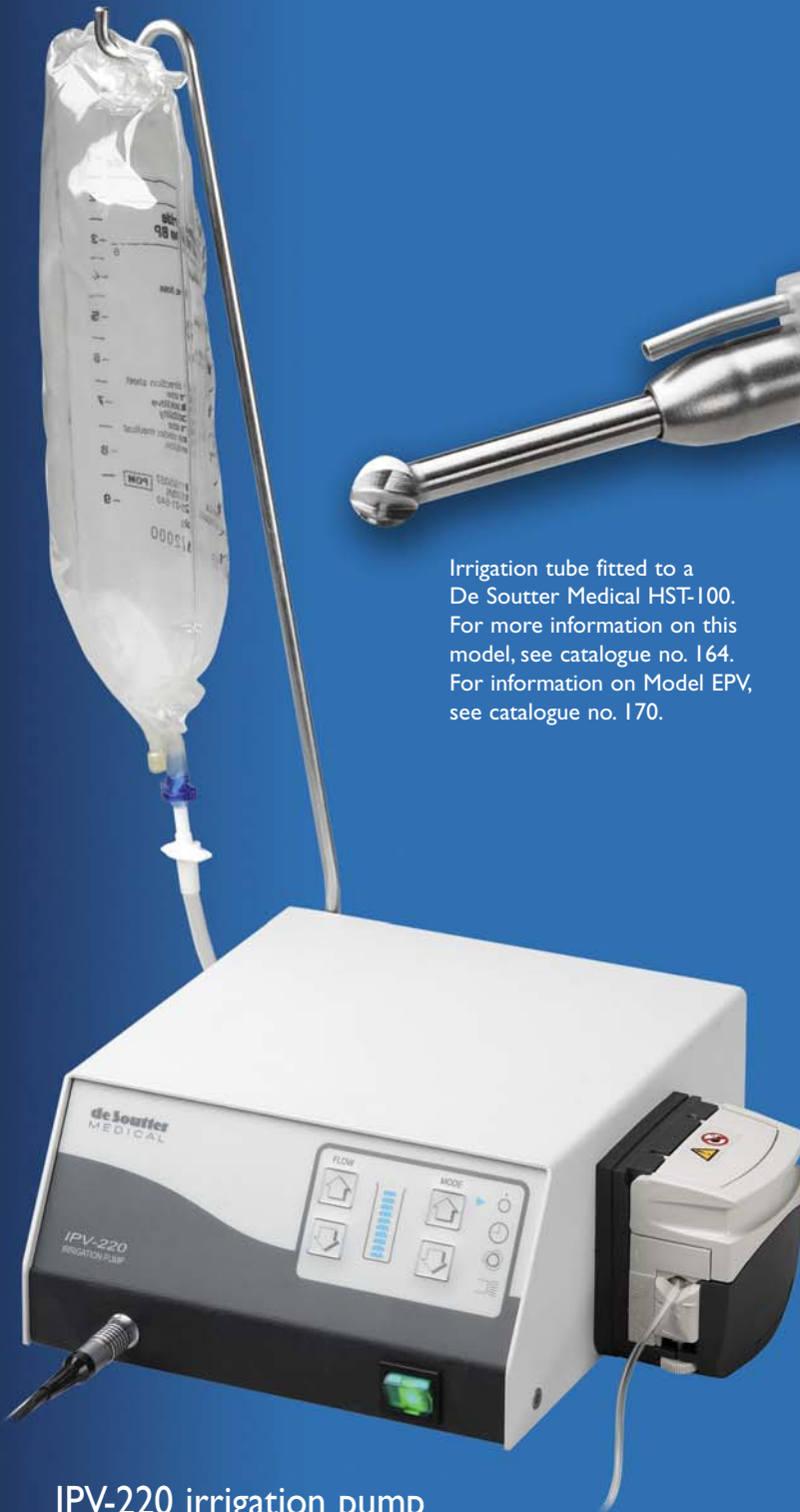

## IPV-220 irrigation pump

100 - 240 V Order code: 11560

Mains power cord sets - 3 m long

UK Order code: 8290

Euro Order code: 8300

US/Canada Order code: 8310

Australia/N Zealand Order code: 8320

Maximum bag capacity - 1 litre. (bag not supplied)

### De Soutter Medical Limited

River Park, Billet Lane, Berkhamsted, Hertfordshire, HP4 1HL. UK

tel: +44 0 1442 860300 fax: +44 0 1442 860333

email: [info@de-souttermedical.com](mailto:info@de-souttermedical.com) web: [www.de-soutter.com](http://www.de-soutter.com)

publication no. 178-1

**de Soutter**  
MEDICAL

# IPV-220 IRRIGATION PUMP

Automatic irrigation system for pneumatic powered instruments used in small bone surgery



## IPS-220 Pump Sensor

Detects small amounts of airflow in the exhaust to start the pump. Fits simply between the air hose and air supply.



## IPV-220 Irrigation Pump

Its design allows the user to quickly and easily select between 4 function settings and adjust flow level. The irrigation tubing clips neatly along the hose to the handpiece to ensure it is kept out of the way.

- Compatible with HST-100 and EPV-220 instruments.
- Fully disposable irrigation kits.
- Simple 'clip-on' set-up.
- 4 function settings.
- Variable flow rate.
- Removeable pump head for cleaning.
- Works with foot or lever control.

IPS-220 pump sensor. Connects to foot pump or air supply via the MA7 probe.

Irrigation tube connected to a De Soutter Medical HST-100 high speed drill.