

dermadrive™

Dermatome skin grafting system


GD SERIES




deSouther
MEDICAL

GD PRECISION SKIN GRAFTING SYSTEM

The slimline and ergonomically designed GD Series dermatome from De Soutter Medical produces consistent, uniform skin grafts time after time. The smooth high speed cutting action ensures an almost effortless and vibration free procedure. The GD-113 dermatome head is designed to be used with the CM2 controller system.

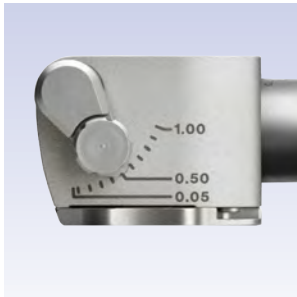

Cut thicknesses are easily selected and constantly maintained using the adjustable depth gauge.

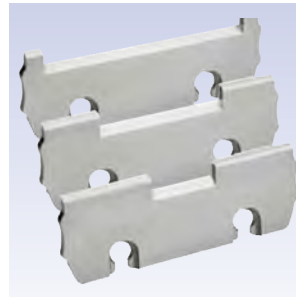

Accurately controlled graft widths of 73 mm, 48 mm and 23 mm are readily achieved by fitting the appropriate easy-fit sole plate.




Stericut precision ground, double facet, ultra sharp blades provide accurate and consistent cutting performance whilst permitting very thin skin grafting. Blades are supplied sterile packed (5 per box). The GD-113 dermatome head enables quick and safe blade replacement.



SYSTEM SPECIFICATION

DERMATOME SPECIFICATION

	Cutting width mm	Operating speed range rpm	Weight kg
	73	2000 - 15000	0.28
The dermatome head is supplied with a 73 mm sole plate as standard. Optional sole plates: 48 mm cutting width - Order code 17260 23 mm cutting width - Order code 17340			

MICRO MOTOR AND POWER SUPPLY

The dermatome head fits quickly onto the universal micro motor which is powered by the variable speed CM2 controller. The motor can be controlled either by the footswitch or directly from the controller itself.

Please specify voltage and plug configuration when ordering.


Power supply specification

Input voltage: 230/115/100 50-60Hz
Power: 60VA
Dimensions: 181 x 101 x 119 mm
Weight: 1.79/0.12 kg controller/motor



STERILISATION CONTAINER

This purpose designed compatible delivery system will ensure your instruments are maintained in good condition at all times. The container is 100% sterilisable by steam.

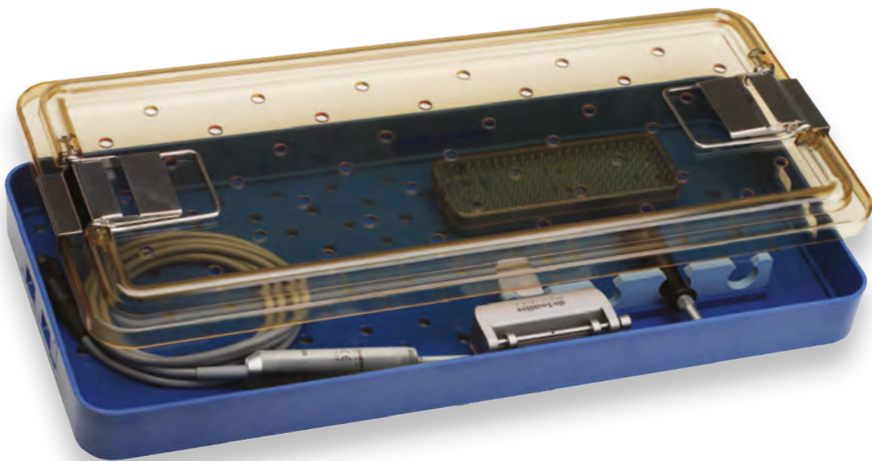

HOW TO ORDER

System components	Order code	Mains power cordsets 3 metres long. Please specify - Order code	
Comprising of contents listed below			
GD-113 head	17200	UK	8290
CM2 motor controller, micro motor and foot pedal - Universal	1292334	Euro	8300
Blades (pack of 5)	S11-103	US/Canada	8310
Sterilisation container	8570	Australia	8320
		New Zealand	8320
		Swiss	14700
Complete system ordering code	GDA-113/230UK		

For further information on De Soutter Medical powered instruments and cast saw systems please contact the address in your area.

UNITED KINGDOM (INTERNATIONAL HEAD OFFICE)

De Soutter Medical Limited
Halton Brook Business Park, Weston Road,
Aston Clinton, Aylesbury, Bucks HP22 5WF
Tel: +44 (0) 1296 634 000
Email: info@de-soutter.com

AUSTRALIA

De Soutter Medical Australia Pty Ltd
2/12-14 Apollo Drive, Hallam
Victoria 3803
Tel: +61 (0) 3 9702 4441
Email: australia@de-soutter.com

AUSTRIA

De Soutter Medical Austria
Hauptstrasse 122-124/Büro 01
3420 Kritzensdorf
Österreich
Tel: +43 (0) 224 32 39000
Email: austria@de-soutter.com

BELGIUM

De Soutter Medical Belgium
Bessemersstraat 14
3620 Lanaken
Tel: +32 (0) 89/47 15 37
Email: belgium@de-soutter.com

FRANCE

De Soutter Medical France
1 252 Avenue du Parc des Expositions
33260 La Teste de Buch
Tel: +33 (0) 5 56 54 89 36
Email: france@de-soutter.com

GERMANY

De Soutter Medical Deutschland
Bahnhofstraße 4
D-66625 Nohfelden
Tel: +49 (0) 68 52-99 12 46
Email: deutschland@de-soutter.com

ITALY

De Soutter Medical Italia
Località Fornace SNC
27022 Casorate Primo - PV
Tel: +39 (0) 2 9009 4098
Email: italy@de-soutter.com

NETHERLANDS

De Soutter Medical Nederland
Gelderlandhaven 2X
3433 PG Nieuwegein
Tel: +31 (0) 85 0491480
Email: nederland@de-soutter.com

UNITED STATES OF AMERICA

De Soutter Medical USA, Inc.
224 Rolling Hill Road, Suite 12A
Mooresville, NC 28117
Tel: +1 (704) 655 9040
Email: usa@de-soutter.com